


PAJUNK®



SonoPlex®

*Echogenic, Stimulating Single Shot
Nerve Block Needles*

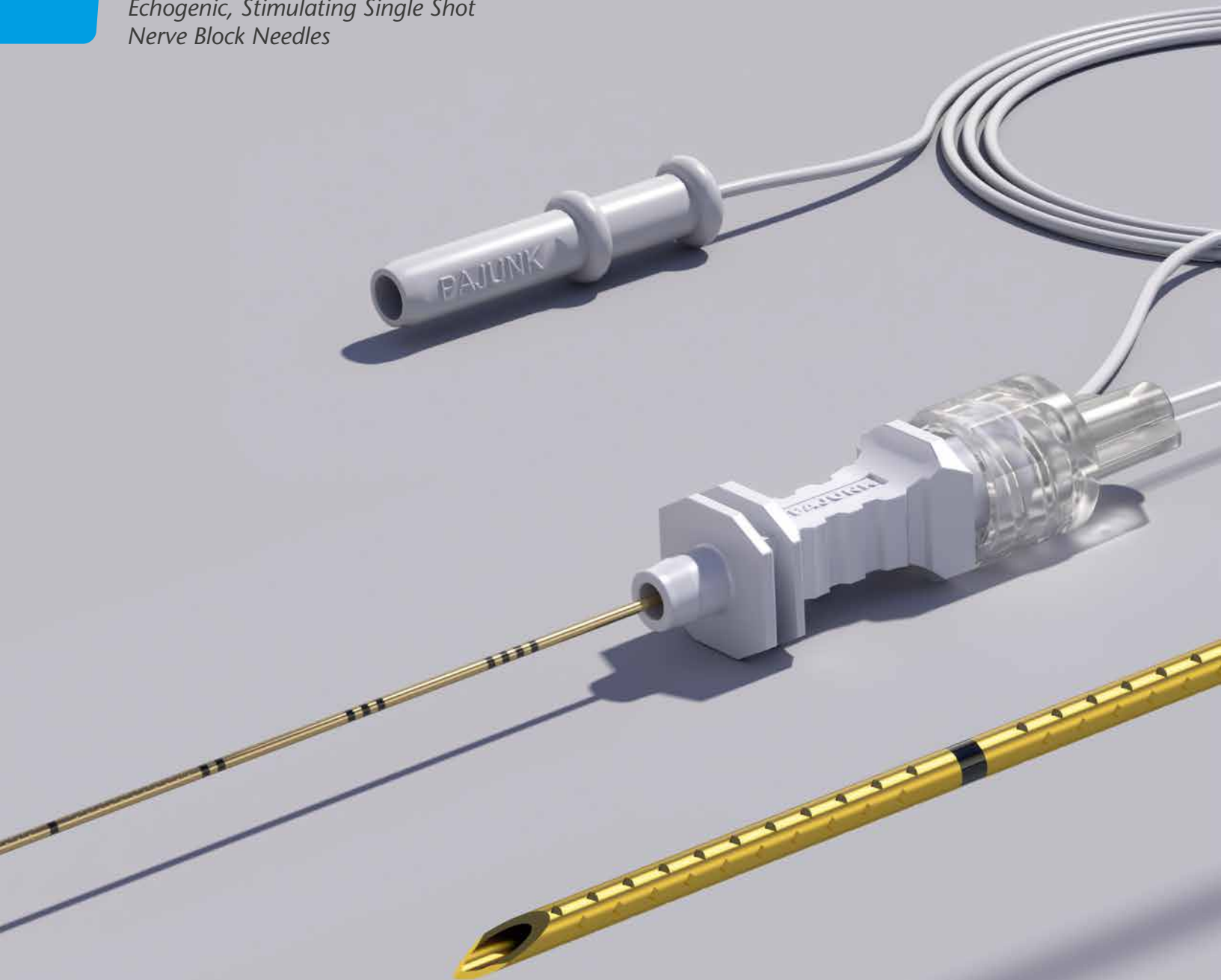

SonoPlex®

The Echogenic Needle with optimized US-Visibility



SonoPlex needles were especially developed by PAJUNK® for single shot nerve block applications, allowing for placement under ultrasound and stimulation techniques (Dual Guidance). As one of the leading echogenic nerve block needles globally, the SonoPlex includes the patented Cornerstone Reflectors for optimized needle echogenicity.



Cornerstone Reflectors

360 degree graduations on the first 20 mm of the needle


- Optimized ultrasound visibility of needle shaft and needle tip¹
- Reliable and optimized needle echogenicity at higher angle

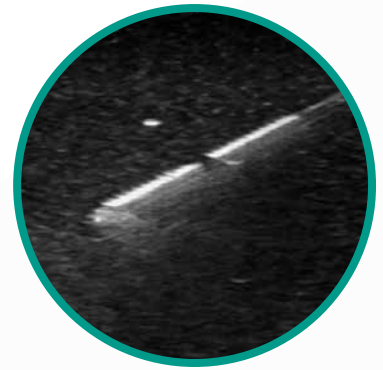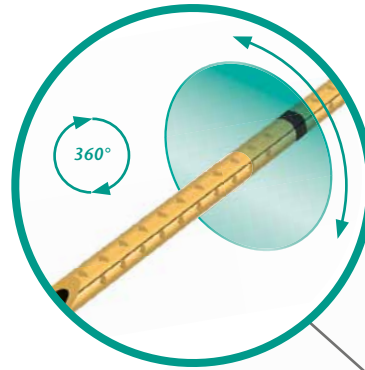

NanoLine® Coating

Very thin polymer film, smooth surface, accurate inner and outer diameter

- Excellent puncture and gliding properties through smooth surface
- Precise stimulation through the non-insulated needle tip

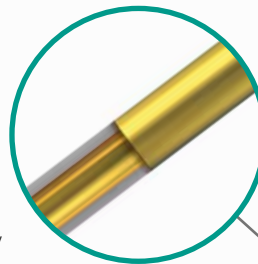

Depth Markings


- Easy to read and identify

Facet tip

Echogenic Needle Tip

Facet tip with two inclination angles

- Improves needle tip visibility under ultrasound

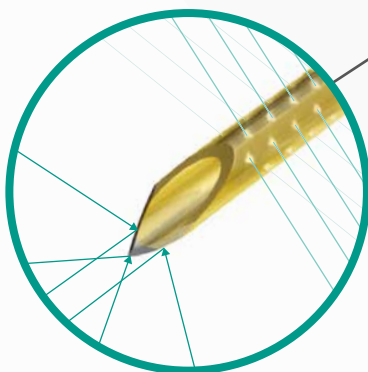

SonoPlex®

SonoPlex needle with facet tip

Size	Item no.	PU
22G x 40 mm (1 5/8")	001185-70	10
22G x 50 mm (2")	001185-74	10
22G x 80 mm (3 1/5")	001185-71	10
21G x 100 mm (4")	001185-77	10
20G x 120 mm (4 3/4")	001185-72	10
20G x 150 mm (6")	001185-76	10

SonoPlex needle with facet S tip

Size	Item no.	PU
22G x 50 mm (2")	001187-74	10
22G x 80 mm (3 1/5")	001187-71	10
21G x 100 mm (4")	001187-77	10
20G x 120 mm (4 3/4")	001187-72	10
20G x 150 mm (6")	001187-76	10

SonoPlex needle with SPROTTE® tip

Size	Item no.	PU
22G x 50 mm (2")	001185-31G	10
22G x 70 mm (3")	001185-31H	10



AVAILABLE OPTIONS

Additional needles for single shot nerve blocks:

- ▶ **SonoBlock II**
Echogenic, non-stimulating single shot nerve block needle, NanoLine® Coating
- ▶ **SonoTAP® II**
Echogenic, non-stimulating single shot plane block needle

More information at: pajunk.com



Scan now!

PAJUNK® GmbH Medizintechnologie
International Sales
Karl-Hall-Straße 1
78187 Geisingen, Germany
Tel.: +49 7704 9291-0
Fax: +49 7704 9291-600
info@pajunk.com
pajunk.com

PAJUNK® Medical Produkte GmbH
Central Europe
Karl-Hall-Straße 1
78187 Geisingen, Germany
Tel.: +49 7704 8008-0
Fax: +49 7704 8008-150
info@pajunk-medipro.com
pajunk.de & pajunk.eu

PAJUNK® Swiss AG
Switzerland
Neuhofstrasse 5A
6340 Baar, Switzerland
Tel.: +41 417 609 627
info@pajunk.ch
pajunk.ch

PAJUNK® Medical System L.P.
United States of America
4575 Marconi Drive
Alpharetta, GA 30005, USA
Tel.: +1 888 972 5865
Fax: +1 678 514 3388
info@pajunk-usa.com
pajunkusa.com

PAJUNK® UK Medical Products Ltd
United Kingdom
Unit D1 The Waterfront
Goldcrest Way, Newburn Riverside
Newcastle upon Tyne, NE15 8NY, UK
Tel.: +44 191 264 7333
info@pajunk.co.uk
pajunk.co.uk